


Understanding the Art of Stenting

EasyRap

Rapidolimus™ Eluting Stent System

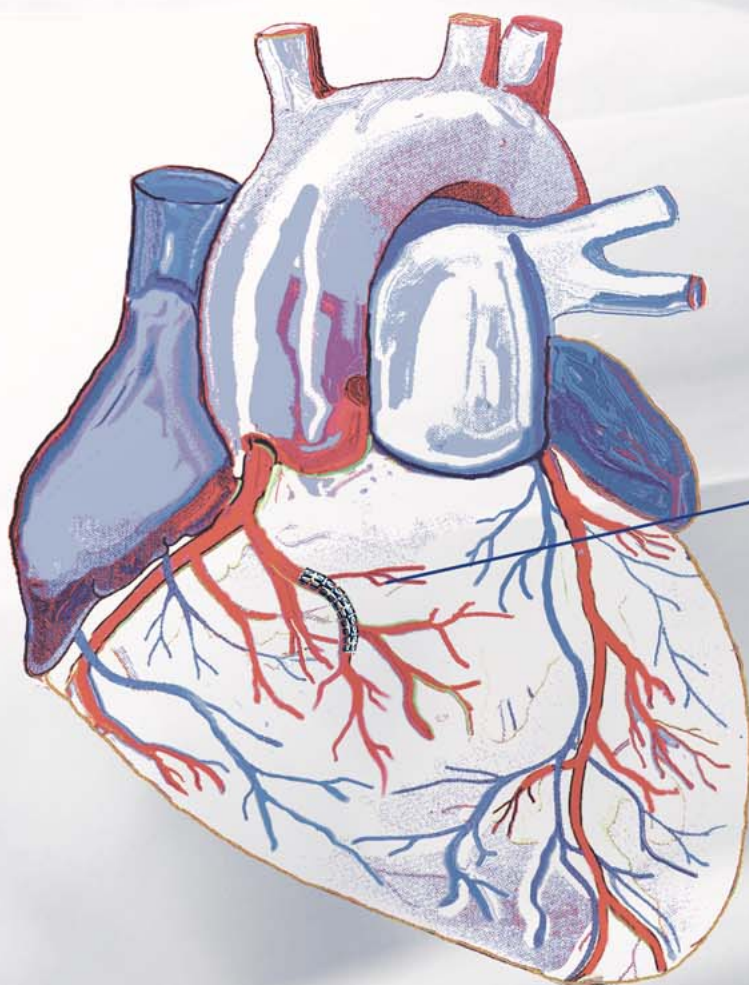

based on the
clinically proven
Easyflex
stent design

EasyRap

Rapidolimus™ Eluting Stent System


EasyRap

Rapidolimus™ Eluting Stent System


Target dose: 1 µg/mm ²
Biosteel coating technology
Material: Stainless steel
DLC (diamond like carbon coated) covered with polymer & drug
Therapeutic release within 28 days
High mechanical stability of the coating

Technical specifications

Drug:	Rapidolimus™ 1µg/mm ²
Polymer:	Inert polymer compound/fully resorbable
Release kinetic:	Slow release within 10-28 days
Coating thickness:	1-2 µm
Stent surface finishing:	DLC PII (plasma ion implantation)

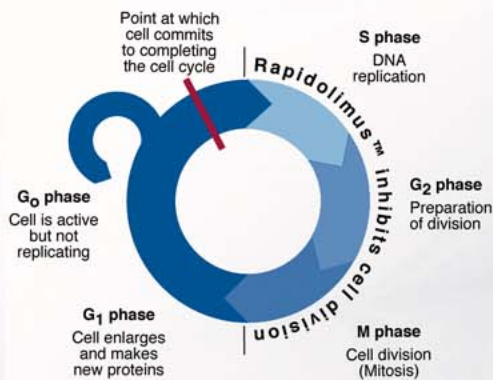


Strut Composition

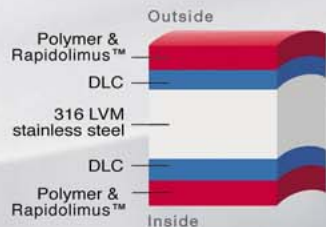

Longitudinal cut through stent strut

Two opposite articulations between every stent segment guarantee multidirectional flexibility

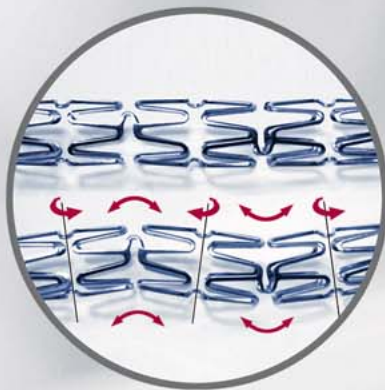

Bare Stent

Catalog Number	Stent Length
ER 0700B	07 mm
ER 1000B	10 mm
ER 1300B	13 mm
ER 1700B	17 mm
ER 2000B	20 mm
ER 2400B	24 mm
ER 2800B	28 mm
ER 3300B	33 mm
ER 3800B	38 mm

Balloon Mounted Stent

Catalog Number	Balloon Diameter	Stent Length
ER 0725	2,50 mm	07 mm
ER 0727	2,75 mm	07 mm
ER 0730	3,00 mm	07 mm
ER 0735	3,50 mm	07 mm
ER 0740	4,00 mm	07 mm
ER 1025	2,50 mm	10 mm
ER 1027	2,75 mm	10 mm
ER 1030	3,00 mm	10 mm
ER 1035	3,50 mm	10 mm
ER 1040	4,00 mm	10 mm
ER 1325	2,50 mm	13 mm
ER 1327	2,75 mm	13 mm
ER 1330	3,00 mm	13 mm
ER 1335	3,50 mm	13 mm
ER 1340	4,00 mm	13 mm
ER 1725	2,50 mm	17 mm
ER 1727	2,75 mm	17 mm
ER 1730	3,00 mm	17 mm
ER 1735	3,50 mm	17 mm
ER 1740	4,00 mm	17 mm
ER 2025	2,50 mm	20 mm
ER 2027	2,75 mm	20 mm
ER 2030	3,00 mm	20 mm
ER 2035	3,50 mm	20 mm
ER 2040	4,00 mm	20 mm
ER 2425	2,50 mm	24 mm
ER 2427	2,75 mm	24 mm
ER 2430	3,00 mm	24 mm
ER 2435	3,50 mm	24 mm
ER 2440	4,00 mm	24 mm
ER 2825	2,50 mm	28 mm
ER 2827	2,75 mm	28 mm
ER 2830	3,00 mm	28 mm
ER 2835	3,50 mm	28 mm
ER 2840	4,00 mm	28 mm
ER 3325	2,50 mm	33 mm
ER 3327	2,75 mm	33 mm
ER 3330	3,00 mm	33 mm
ER 3335	3,50 mm	33 mm
ER 3340	4,00 mm	33 mm
ER 3825	2,50 mm	38 mm
ER 3827	2,75 mm	38 mm
ER 3830	3,00 mm	38 mm
ER 3835	3,50 mm	38 mm
ER 3840	4,00 mm	38 mm

distributed by

iramed medical devices

iramed GmbH, Germany
 Im End 39, D-41844 Wegberg-Arsbeck
 Tel. +49-2436-588, Fax +49-2436-1091
 E-mail: iramed@aol.com, www.iramed.com